

We'll lead in **innovation**
We'll lead in **efficiency**
We'll lead in **safety**

BestoCure
For their well-being and future

wellead | BestoCure

WELL LEAD MEDICAL CO., LTD.

Add: 47 Guomao Avenue South, Panyu, Guangzhou, China

Tel: +86-20-84758878 Fax: +86-20-84758224

Post Code: 511434

E-mail: info@wellead.com.cn

www.wellead.com

Veterinary Catalogue

For their well-being and future.

We are more than just a manufacturer. We bring the same values and standards from our global spreading human healthcare solutions into the veterinary field—offering smarter, more compassionate care for both animals and the professionals who treat them.

As an established manufacturer, we partner with leading medical device companies worldwide. We are committed to delivering products that meet the high standards of quality. With 27 years of supplying expertise, an R&D team of over 400 employees, and an inner national-level laboratory, we are also here to listen and bring your OEM product ideas to life.

Our Vision

Strive to advance the animal health care by being the world-wide leading provider of veterinary solutions and services.

wellead | **BestoCure**

www.wellead.com



PRODUCT
SOLUTIONS

Urology

- Cat Catheterisation
- Silicone Foley Catheters
- Dog Catheterisation
- Urine Collection
- Accessories

Anesthesia

- Airway Management
- Pain Management
- Regional Anesthesia

Hemodialysis

- Blood Circuits
- Hemodialysis Catheters

Respiratory

- Aerosol Chamber
- Oxygen Cannula

Wound Management

- Drainage System & Tubes

Nutrition

- Feeding Tubes

Vascular Access


- Central Venous Catheters

AKVOGLIDE™

Hydrophilic Series

Hydrophilic coated catheters have a layer of a pre-lubricated coating that is bound to the catheter surface. The pre-lubricated coating absorbs and binds water to the hydrophilic catheter, resulting in a thick, smooth, and slippery surface. The coating layer remains intact upon introduction into the urethra and ensures lubrication of the urethra in its entire length.

Experience the glide.
Redefining smooth.

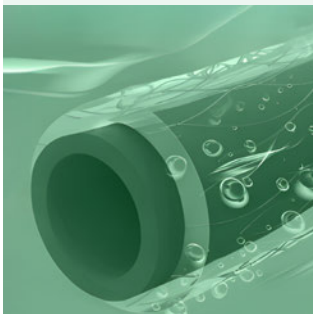

MAIN BENEFITS AT A GLANCE

- ① Effortless catheter insertion with 99% less friction than conventional gel-based solutions.
- ② Maintains lubrication throughout indwelling, eliminating the risk of drying.
- ③ Smooth insertion and removal minimize discomfort and irritation, particularly in tomcat urethra.
- ④ Instant activation with water.
- ⑤ Extensively validated through clinical use in human healthcare.

Performance of Hydrophilic Coating		
Performance indicators	Wellead products	Clinical significance
Friction coefficient	≤0.05 (Gel-lubricated: >0.5)	Patients experience comfortable catheterisation, reducing the occurrence of difficult catheterisations
Adhesion firmness & durability	After 30 days of immersion in artificial urine and 30 repeated friction cycles, the coefficient of friction decreased by 5%	The catheter was removed smoothly after 14 days of indwelling
Activation time	<5 seconds	Reduces healthcare professionals' procedural time
Surface drying time	≥8 mins	Eliminates the need for repeated lubrication

Trusted worldwide in human healthcare, Wellead's exclusive-patent AkvoGlide™ hydrophilic coating sets the standard with its exceptionally low friction, superior durability, and flawless uniformity—engineered for unparalleled performance.

Proven in veterinary practice, AkvoGlide™ directly improves patient outcomes. It significantly lowers the risk of UTIs and microscopic hematuria in catheterized cats, offering gentle protection that minimizes irritation and mechanical damage to the delicate urethral and bladder lining throughout indwelling.


Start with us to share the same welfare as human healthcare.

AKVOGLIDE™

Hydrophilic Cat Urinary Catheterisation Set

Your new catheterisation solution: catheters for blockage flushing and indwelling both in one integrated pack.

- **Components**
 - AkvoGlide™ Cat Urinary Catheter, side holes, with stylet
 - Flushing Catheter, open ended
 - Sterile water sachet
 - Lubrication gel
 - Surgical drape
 - Inner beveled tray

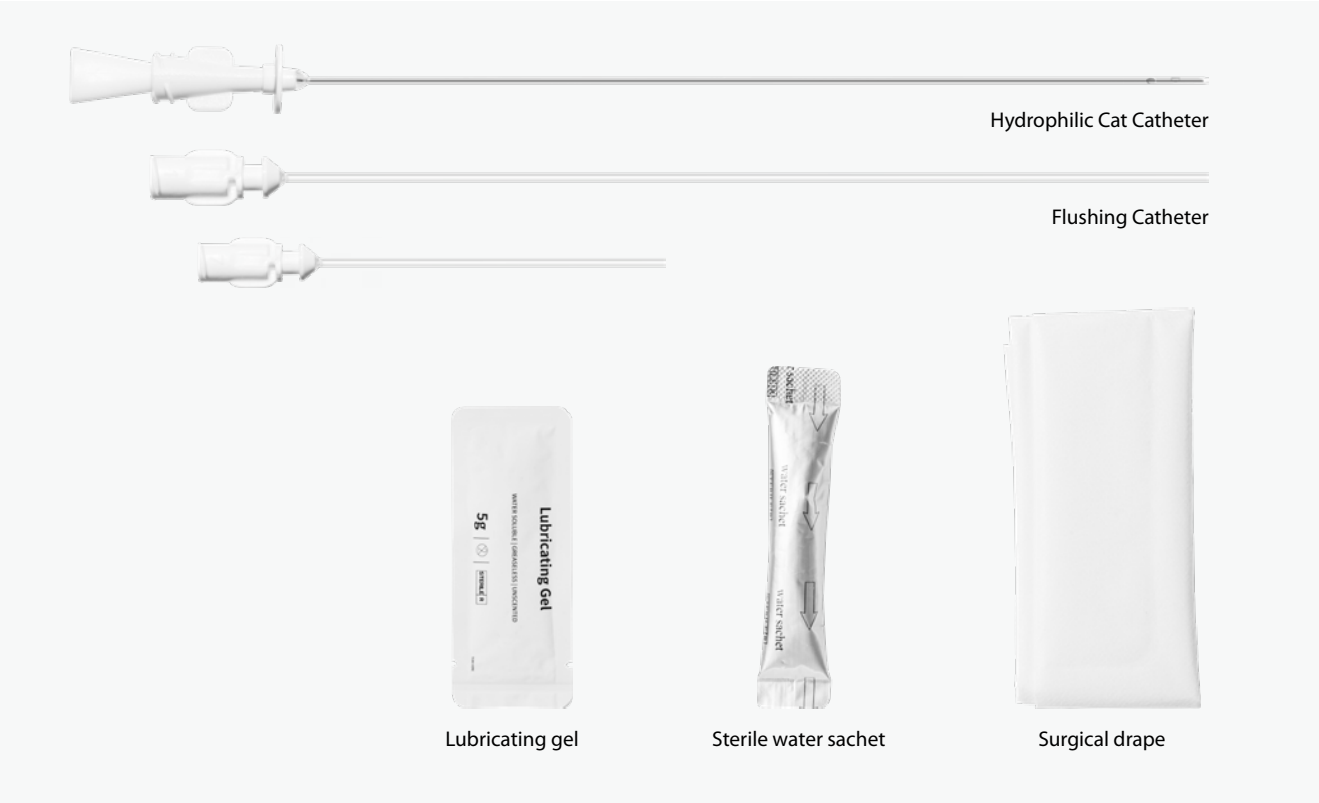


Code	Indwelling Catheter OD X Length	Flushing Catheter OD X Length	Pack
VU04S0310A	3Fr 1.0 x 100mm	3Fr 1.0 x 50mm	10 pcs/box
VU04S0313A	3Fr 1.0 x 130mm	3Fr 1.0 x 50mm	
VU04S0413A	4Fr 1.3 x 130mm	3Fr 1.0 x 50mm	
VU04S0310B	3Fr 1.0 x 130mm	3Fr 1.0 x 130mm	
VU04S0413B	4Fr 1.3 x 130mm	3Fr 1.0 x 130mm	

Hydrophilic Cat Urinary Catheter

Side holes, with stylet

- Easy to activate by simply adding sterile water or saline in the hard blister package groove after peeling up
- With superior biocompatibility, long term placement(up to 10 days)
- Ultra-soft with effective resistance to kinking, improving comfort degree
- With stainless steel stylet, ensuring easy catheterisation

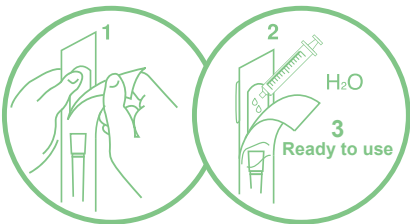

HOW TO USE




Code	Specifications	Pack
VU040310	3Fr 1.0 x 100mm	12 pcs/box
VU040313	3Fr 1.0 x 130mm	
VU040413	4Fr 1.3 x 130mm	
VU040416	4Fr 1.3 x 160mm	

Hydrophilic Silicone Foley Catheter

2-way, transparent, with stylet

- **Accessories**
- Foley-Luer adaptor
- Sterile water sachets
- Tray

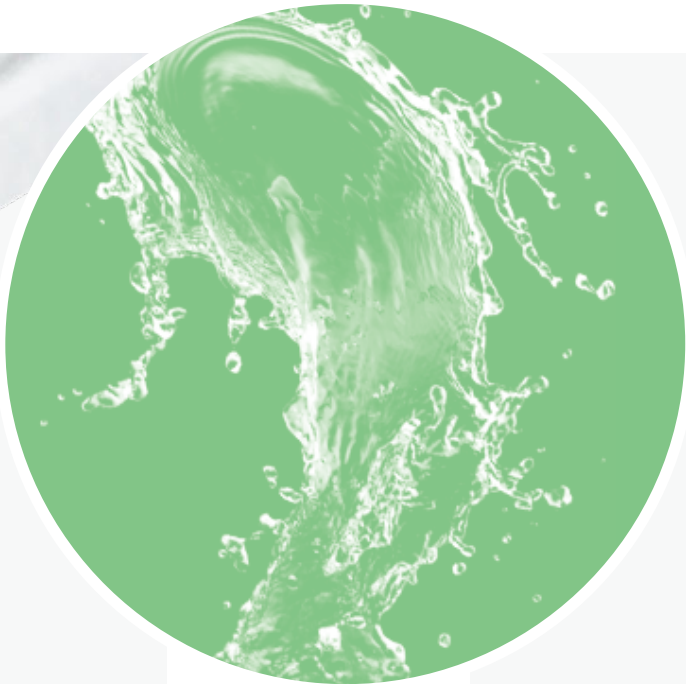

Code	Balloon Size (ml)	Size (Fr)	Length (cm)	Pack
VU01H0630	3	6	30	5 pcs/box
VU01H0655	3	6	55	
VU01H0670	3	6	70	
VU01H0830	3-5	8	30	
VU01H0855	3-5	8	55	
VU01H0870	3-5	8	70	
VU01H1055	3-5	10	55	
VU01H1070	3-5	10	70	

AKVOGLIDE™

Hydrophilic Dog Catheter

Side holes

- Accessories: sterile water sachet

Code	Specifications	Pack
VU070650	6Fr 2.0 x 500mm	5 pcs/box
VU070850	8Fr 2.6 x 500mm	
VU071050	10Fr 3.3 x 500mm	



AKVOGLIDE™

Hydrophilic Nasogastric Feeding Tube

- Accessories: suture wing, sterile water sachet

Code	Specifications	Pack
VG02B0650	6Fr x 500mm	10 pcs/box
VG02B0850	8Fr x 500mm	
VG02B1050	10Fr x 500mm	



HOW TO USE

Cat Urinary Catheter (Soft Type)

Side holes, soft suture collar, with stylet

- Designed for less trauma and infection indwelling use
- With superior biocompatibility, long term placement(up to 10 days)
- Ultra-soft with effective resistance to kinking. , improving comfort degree
- With stainless steel stylet, ensuring easy catheterization
- Soft suture collar for higher indwelling comfort degree

Code	Specifications	Pack
VU030310A	3Fr 1.0 x 100mm	12 pcs/box
VU030313A	3Fr 1.0 x 130mm	
VU030413A	4Fr 1.3 x 130mm	
VU030416A	4Fr 1.3 x 160mm	

Cat Urinary Catheter (Soft Type)

Side holes, rigid suture collar, with stylet

- Designed for less trauma and infection indwelling use
- With superior biocompatibility, long term placement(up to 10 days)
- Ultra-soft with effective resistance to kinking. , improving comfort degree
- With stainless steel stylet, ensuring easy catheterization

Code	Specifications	Pack
VU030310	3Fr 1.0 x 100mm	12 pcs/box
VU030313	3Fr 1.0 x 130mm	
VU030413	4Fr 1.3 x 130mm	
VU030416	4Fr 1.3 x 160mm	

Cat Urinary Catheter (Rigid Type)

Side Holes, open end, without stylet

- Designed for better flushing the urinary track and release blockage
- With higher rigidity
- Atraumatic open end

Code	Specifications	Pack
VU050310	3Fr 1.0 x 100mm	20 pcs/box
VU050313	3Fr 1.0 x 130mm	
VU050413	4Fr 1.3 x 130mm	
VU050416	4Fr 1.3 x 160mm	



All Silicone Foley Catheter

2-way, transparent, with stylet

- With great biocompatibility, long term placement for up to 29 days
- With radiopaque line
- With easily removable stainless steel wire stylet (6/8/10 Fr), ensuring easy catheterization

Code	Balloon Size (ml)	Size (Fr)	Length (cm)	Pack
VU010630	3	6	30	5 pcs/box
VU010655	3	6	55	
VU010670	3	6	70	
VU010830	3-5	8	30	
VU010855	3-5	8	55	
VU010870	3-5	8	70	
VU011055	3-5	10	55	
VU011070	3-5	10	70	
VU011255	10	12	55	
VU011455	10	14	55	
VU011655	10	16	55	
VU011855	30	18	55	
VU012065	30	20	65	
VU012265	30	22	65	




Premium All Silicone Foley Catheter

2-way, white, without stylet

- With great biocompatibility, long term placement for up to 29 days
- Whole tubing distinct X-ray detectable
- Increasing rigidity for insertion without wire wire

Code	Balloon Size (ml)	Size (Fr)	Length (cm)	Pack
VU020630	3	6	30	5 pcs/box
VU020655	3	6	55	
VU020670	3	6	70	
VU020830	3-5	8	30	
VU020855	3-5	8	55	
VU020870	3-5	8	70	
VU021055	3-5	10	55	
VU021070	3-5	10	70	
VU021255	10	12	55	
VU021455	10	14	55	
VU021655	10	16	55	
VU021855	30	18	55	
VU022065	30	20	65	
VU022265	30	22	65	

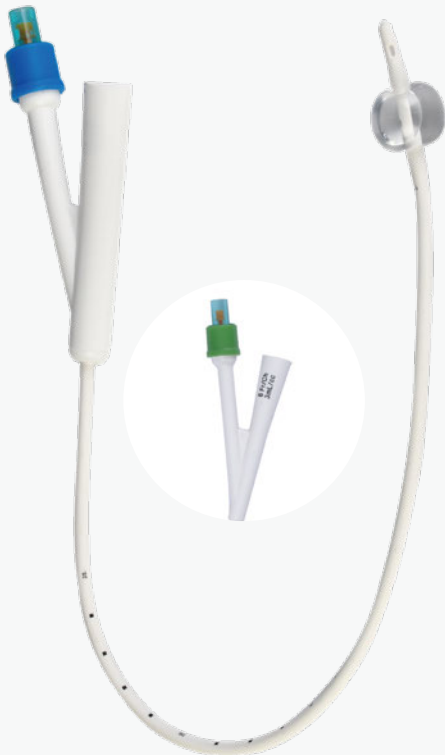


Dog Catheter

Side Holes, without stylet

- Designed for single-time urine drainage or bladder lavage
- Altraumatic closed end
- Female Luer lock fitting
- Colour-coded for size identification

Code	Specifications	Pack
VU060650	6Fr 2.0 x 500mm	10 pcs/box
VU060850	8Fr 2.6 x 500mm	
VU061050	10Fr 3.3 x 500mm	

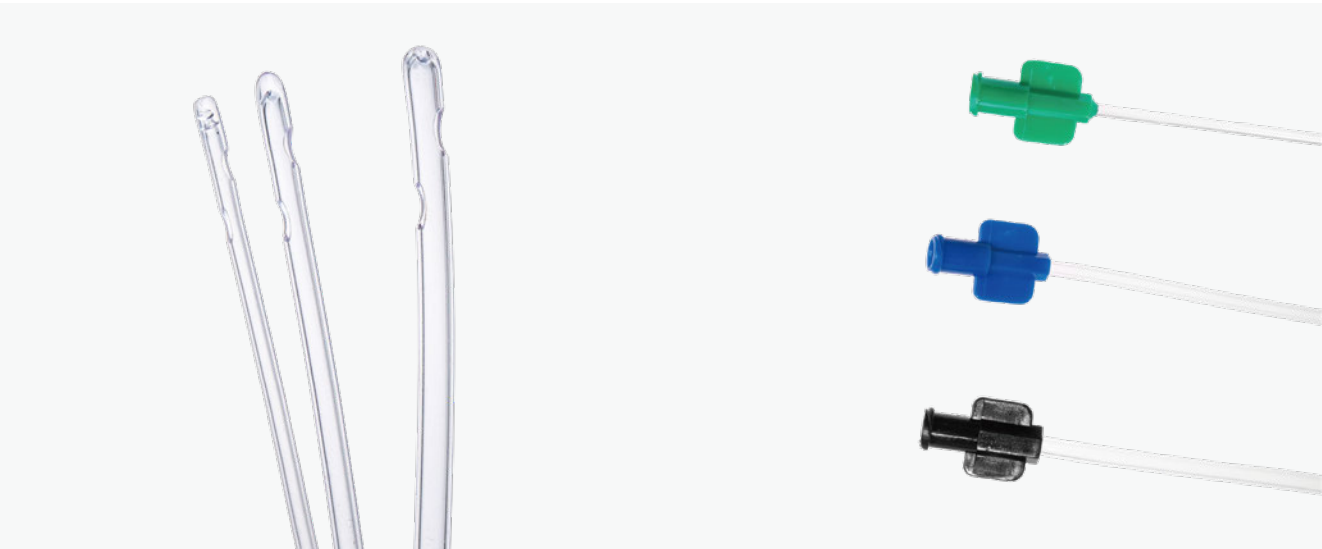

Urine Bag

200/300/500ml

- Non-return valve prevents risk of infection
- With 2 comfort straps. Designed for both pets portable and cage fixing use.

Code (T-Tap Valve)	Code (Flip Valve)	Volume	Connector	Extension Line	Pack
VU08A101	VU08B101	200ml	Luer lock	1m, straight	10 pcs/ bag
VU08A111	VU08B111		Luer lock	1.5m, sprial	
VU08A201	VU08B201		Foley hub	-	
VU08A103	VU08B103	300ml	Luer lock	1m, straight	
VU08A113	VU08B113		Luer lock	1.5m, sprial	
VU08A203	VU08B203		Foley hub	-	
VU08A102	VU08B102	500ml	Luer lock	1m, straight	
VU08A112	VU08B112		Luer lock	1.5m, sprial	
VU08A202	VU08B202		Foley hub	-	



Urinary Catheterization Kit

- Containing different urinary catheters and other components for catheter insertion or exchange
- Components can be customized to meet local guidelines and different using habits

- **Components**
- Standard
- Urinary Catheter, Urine Bag, Medical Glove



- **Optional**
- Spigot

Foley-Luer Adaptor

Stasis Valve

Empty Syringe

Sterile Water Prefilled Syringe

Lubricant Gel Prefilled Syringe
- Lubricant Gel Bag

3-0 Nylon Suture

Suture Wing with clamp

Wrap

Underpad

Fenestrated Drape
- Gauze

Medical Absorbent Cotton

Environment-friendly Tray

Povidon Iodine Swabsticks

Alcohol Pad

*Available Catheters: Cat Catheters, Hydrophilic Cat Catheters, All Silicone Foley Catheters, Hydrophilic Silicone Foley Catheters, All Silicone Premium Foley Catheters

*Available Urine bags: 200ml/500ml/750ml/1000ml, Luer/Foley Fitting

Accessories

Spigot Catheter(Catheter Plug)



Foley-Luer Adaptor(Catheter Cap)

To provide stop of flow for catheters hygienically or transform the Foley drainage lumen to connect to a urine bag or syringe with Luer lock fitting during nursing procedures.

- Universal size to fit all sizes of Foley Catheters
- Ribbed design provides strong and secure connection

Code	Description
VU09A0101	Spigot Catheter
VU09A0102	Foley-Luer Adaptor




Lubricant Gel

- Packet (5g)
- Pre-filled syringe

To provide efficient lubricating support during catheters inserting procedures.

- 100% Water-based. Water soluble, greaseless, unscented and non-irritating gel for various surgical lubrication.

Code	Description
VU10A0105	Packet (5g)
VU10B0105	Pre-filled syringe (5ml)

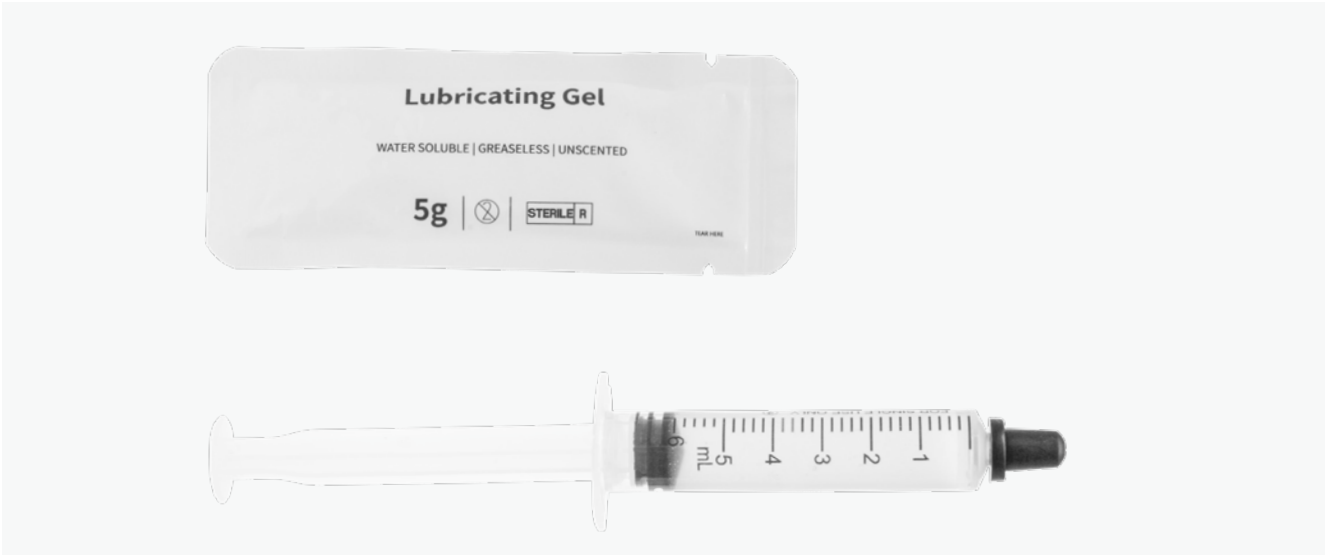

Autoclavable Endotracheal Tube, Silicone

- Transparent, soft, with radiopaque line
- Made of medical grade silicone with great biocompatibility
- Reusable, Can be sterilized in the autoclave at 121°C (250°F) up to 20 times

Code Standard Tip	Size ID(mm)	Pack
VA02A030	3.0	1 pc/pk
VA02A040	4.0	
VA02A050	5.0	
VA02A060	6.0	
VA02A070	7.0	
VA02A080	8.0	
VA02A090	9.0	
VA02A100	10.0	
VA02A110	11.0	
VA02A120	12.0	
VA02A140	14.0	



Endotracheal Tube, Silicone, Avian/Exotic Animals

- Both with flexibility and rigidity
- Larger proximal for tighter seal
- With guide wire (ID 1.0mm & 1.5mm) , providing support to intubation

Code Standard Tip	Size ID(mm)	Pack
VA03A010	1.0	1 pc/pk
VA03A015	1.5	
VA03A020	2.0	
VA03A025	2.5	
VA03A030	3.0	
VA03A035	3.5	
VA03A040	4.0	
VA03A050	5.0	
VA03A060	6.0	



Veterinary Autoclavable Endotracheal Tube, Silicone, Foal

Code Standard Tip	Size ID(mm)	Length(mm)	Pack
VA04A070	7.0	550	1 pc/pk
VA04A080	8.0	550	
VA04A090	9.0	550	
VA04A100	10.0	550	
VA04A110	11.0	550	
VA04A120	12.0	550	
VA04A140	14.0	550	

Standard Endotracheal Tube, PVC

Side Holes, without stylet

- **Two tip types:**
 - Standard tip
 - Bull nose tip, for less trauma intubation
- Available with DEHP-Free material
- Available with pre-loaded intubating stylet
- Transparent, kink-resistant, with radiopaque line
 - Murphy Eye: reducing the risk of occlusion and maintaining airflow

Uncuffed

Code Standard Tip	Code Bull Nose Tip	Size ID (mm)	Pack
VA01A120	VA01A220	2.0	10 pcs/box
VA01A125	VA01A225	2.5	
VA01A130	VA01A230	3.0	
VA01A135	VA01A235	3.5	
VA01A140	VA01A240	4.0	
VA01A145	VA01A245	4.5	
VA01A150	VA01A250	5.0	
VA01A155	VA01A255	5.5	
VA01A160	VA01A260	6.0	
VA01A165	VA01A265	6.5	
VA01A170	VA01A270	7.0	

High Volume Low Pressure Cuff (Columnar Cuff)

Code Standard Tip	Code Bull Nose Tip	Size ID (mm)	Pack
VA01B120	VA01B220	2.0	10 pcs/box
VA01B125	VA01B225	2.5	
VA01B130	VA01B230	3.0	
VA01B135	VA01B235	3.5	
VA01B140	VA01B240	4.0	
VA01B145	VA01B245	4.5	
VA01B150	VA01B250	5.0	
VA01B155	VA01B255	5.5	
VA01B160	VA01B260	6.0	
VA01B165	VA01B265	6.5	
VA01B170	VA01B270	7.0	
VA01B175	VA01B275	7.5	
VA01B180	VA01B280	8.0	
VA01B185	VA01B285	8.5	
VA01B190	VA01B290	9.0	
VA01B195	VA01B295	9.5	
VA01B196	VA01B296	10.0	
VA01B197	VA01B297	10.5	
VA01B198	VA01B298	11.0	



Low Profile Cuff (Olive-shaped Cuff)

Code Standard Tip	Code Bull Nose Tip	Size ID (mm)	Pack
VA01C125	VA01C225	2.5	10 pcs/box
VA01C130	VA01C230	3.0	
VA01C135	VA01C235	3.5	
VA01C140	VA01C240	4.0	
VA01C145	VA01C245	4.5	
VA01C150	VA01C250	5.0	
VA01C155	VA01C255	5.5	
VA01C160	VA01C260	6.0	
VA01C165	VA01C265	6.5	
VA01C170	VA01C270	7.0	
VA01C175	VA01C275	7.5	
VA01C180	VA01C280	8.0	
VA01C185	VA01C285	8.5	
VA01C190	VA01C290	9.0	
VA01C195	VA01C295	9.5	
VA01C196	VA01C296	10.0	



Endotracheal Tube, Cuffed, PVC, Larger Sizes


- Curved tubing for easier intubation

Code Standard Tip	Size ID(mm)	Pack
VA01E120	12.0	5 pcs/box
VA01E130	13.0	
VA01E140	14.0	

PU Ultra-thin Cuff (Columnar Cuff)

- Made from PU in Medical grade.
- Reducing intubation-related pneumonia.
- High degree of comfort and soft.Protecting more efficiently against fluid leakage or microaspiration.


Code Standard Tip	Code Bull Nose Tip	Size ID (mm)	Pack
VA01D130	VA01D230	3.0	10 pcs/box
VA01D135	VA01D235	3.5	
VA01D140	VA01D240	4.0	
VA01D145	VA01D245	4.5	
VA01D150	VA01D250	5.0	
VA01D155	VA01D255	5.5	
VA01D160	VA01D260	6.0	
VA01D165	VA01D265	6.5	
VA01D170	VA01D270	7.0	
VA01D175	VA01D275	7.5	
VA01D180	VA01D280	8.0	
VA01D185	VA01D285	8.5	
VA01D190	VA01D290	9.0	
VA01D195	VA01D295	9.5	
VA01D196	VA01D296	10.0	

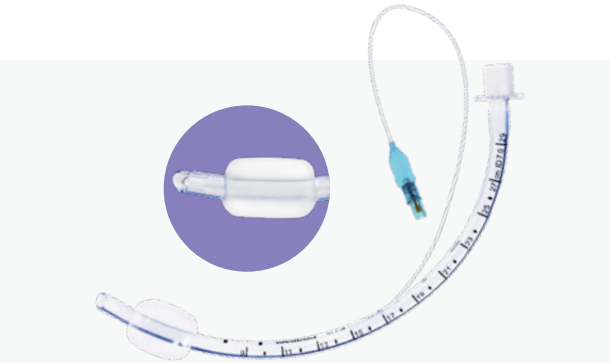

Reinforced Endotracheal Tube, Cuffed, PVC

- With a spiral wire embedded in the tracheal tube provides effective resistance to kinking or bite.
- Better airway for dental works, eye surgery procedures and neurosurgical procedures, etc.

Code Standard Tip	Code Bull Nose Tip	Size ID (mm)	Pack
VA05B130	VA05B230	3.0	5 pcs/box
VA05B135	VA05B235	3.5	
VA05B140	VA05B240	4.0	
VA05B145	VA05B245	4.5	
VA05B150	VA05B250	5.0	
VA05B155	VA05B255	5.5	
VA05B160	VA05B260	6.0	
VA05B165	VA05B265	6.5	
VA05B170	VA05B270	7.0	
VA05B175	VA05B275	7.5	
VA05B180	VA05B280	8.0	
VA05B185	VA05B285	8.5	
VA05B190	VA05B290	9.0	
VA05B195	VA05B295	9.5	
VA05B196	VA05B296	10.0	

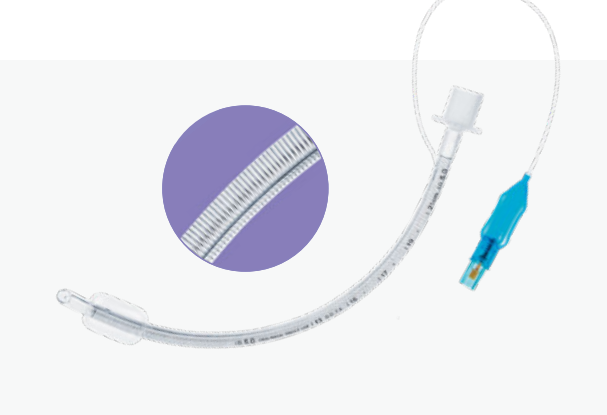

Intubating Stylet

- Made of malleable aluminum covered with plastic sheath.
- Soft rounded tip for minimizing patient trauma during intubation.
- Flexible use for assistance of intubation.

Code	Size	Pack
VA07A06F	6Fr	10 pcs/box
VA07A10F	10Fr	
VA07A12F	12Fr	
VA07A14F	14Fr	

SuctionPlus Endotracheal Tube with Evacuation Lumen, Cuffed, PVC

- SuctionPlus tracheal tube is designed with the evacuation
- lumen in order to provide protection against risk of aspiration and decrease the rate of ventilation-associated pneumonia (VAP).


Code Standard Tip	Code Bull Nose Tip	Size ID (mm)	Pack
VA06B160	VA06B260	6.0	5 pcs/box
VA06B165	VA06B265	6.5	
VA06B170	VA06B270	7.0	
VA06B175	VA06B275	7.5	
VA06B180	VA06B280	8.0	
VA06B185	VA06B285	8.5	
VA06B190	VA06B290	9.0	



Low Dead Space Endotracheal Tube Connector

- Designed for the accurate measurement of end-tidal CO₂ in low-body-weight animals.

Code	Size ID (mm)	Pack
VA09A020	2.0	1 set/pk
VA09A025	2.5	
VA09A030	3.0	
VA09A035	3.5	
VA09A040	4.0	
VA09A045	4.5	

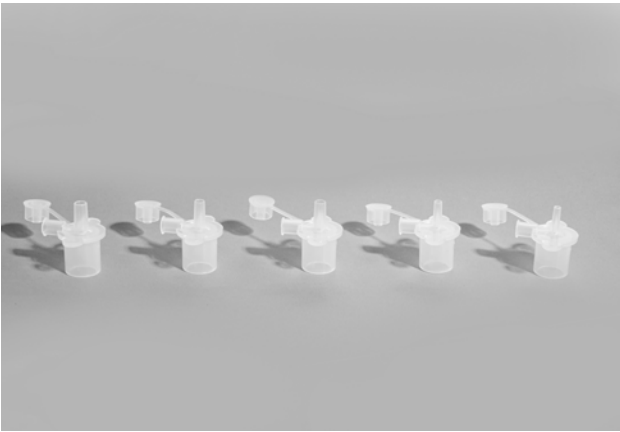

Anesthetic Mask

- Made of clear plastic, with black rubber diaphragm.
- With elastic neck strap (except for #1 & #2) for better fixation of the mask. 15 mm connection piece fits inner patient hose.

Code	Size	ID (mm)	OD (mm)	Pack
VA08A001	1	27	51	Without strp, 15mm connector
VA08A002	2	32	61	
VA08A003	3	34	66	With strap,15mm connector
VA08A004	4	32	93	
VA08A005	5	44	113	With strap and 22mm conector-adaptor, 15mm connector
VA08A006	6	73	138	



Tracheostomy Tube

- Designed for effectively inserting directly into the trachea through a surgical incision in case of emergency first-aid.
- Straps provided in the pack for fixation of the tube. Smooth rounded tip of the Obturator reduces trauma during insertion.

Code	Size ID (mm)	Pack
VA10A030	3	3 pcs/box
VA10A050	5	
VA10A070	7	




Endobronchial Blocker Tube

Endobronchial blocker tube is a device which can be inserted down a tracheal tube after tracheal intubation so as to block off the right or left main bronchus of the lungs in order to be able to achieve a controlled one lung ventilation. The lung tissue distal to the obstruction will collapse, thus allowing the surgeon's view and access to relevant structures within the thoracic cavity.

- Bronchoscope port: Soft and elastic material for safe insertion of a bronchoscope. Twist lock design: Easy operation & secure locking.
- PU Ultra-thin cuff: Creating a superior seal and reduces microaspiration. When deflated, it is very close to the blocker tube, so it can pass ETT easily.

Code	Size	Pack
VA11A05F	5Fr	1 pcs/box
VA11A07F	7Fr	
VA11A09F	9Fr	

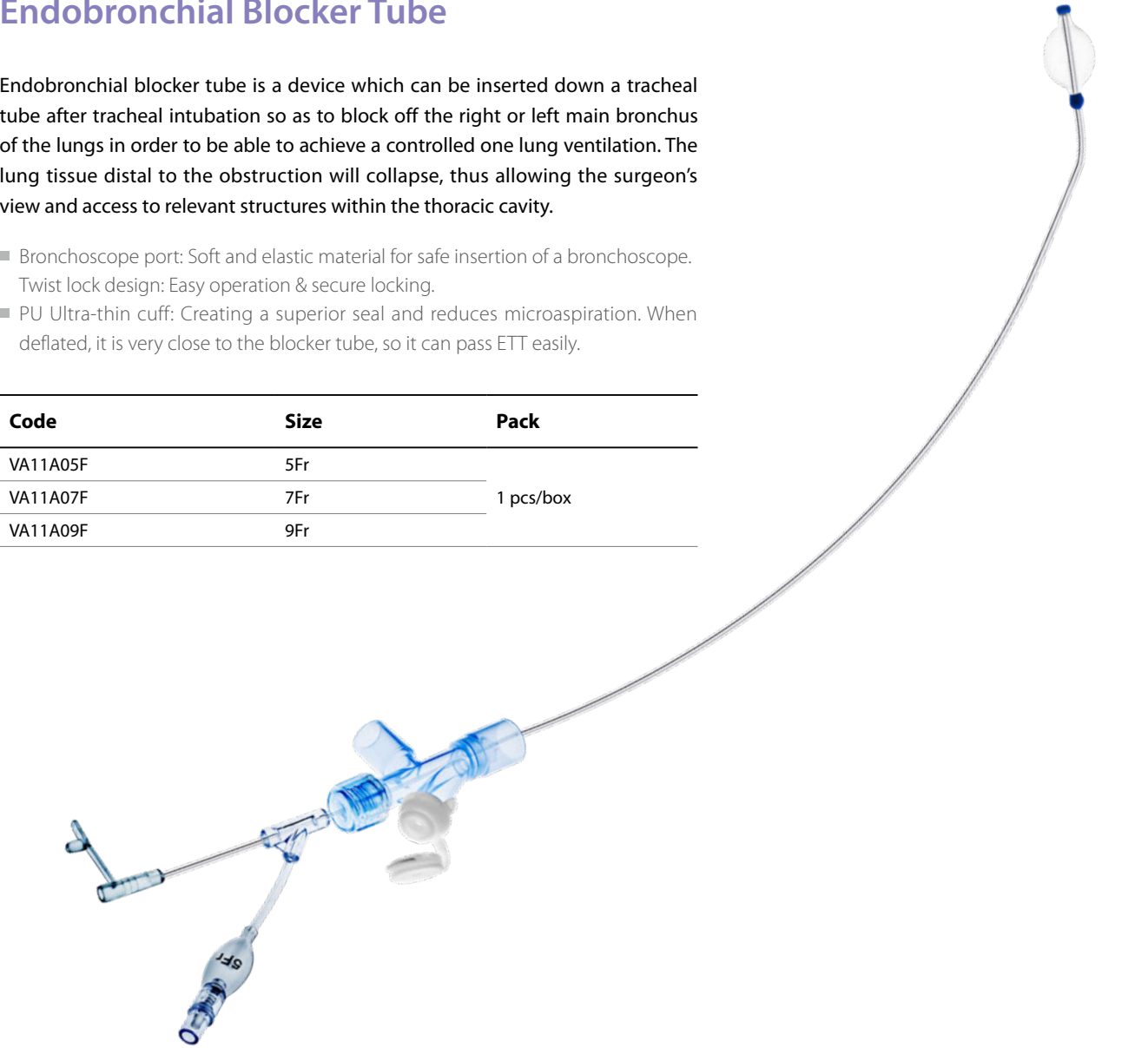

Anesthesia Breathing Circuit, Rebreathing

- Light weight reduces torque on endotracheal tube, etc.
- High transparency. High resistance to collapse
- Corrugated, extendable and Bilumen tubes for choosing
- Length: 1.0m/1.5m/1.8m/2.4m/3.6m



Divisiθ® Bilumen Circuit



Extendable Circuit



Corrugated Circuit

Veterinary Anaesthesia Breathing Circuit, Non-Rebreathing

- Tail-end control valve allows the adjustment of pressure levels within the system
- *With two types:*
 - Undetachable 0.5L breathing bag
 - Detachable 0.25L breathing bag



Breathing Bags, Latex-Free

- Light weight design, ISO standard connectors

Code	Volume(L)	Pack
VA14A005	0.5	10 pcs/box
VA14A010	1	
VA14A020	2	
VA14A030	3	



Breathing Bags, Latex

- Light weight design, ISO standard connectors

Code	Volume(L)	Pack
VA14B050	5	10 pcs/box
VA14B080	8	
VA14B100	10	
VA14B120	12	
VA14B150	15	
VA14B300	30	



Code	Volume(L)	Pack
VA15A012	1.2m, Male - Male Connectors	50 pcs/box
VA15A015	1.5m, Male - Male Connectors	
VA15A020	2.0m, Male - Male Connectors	
VA15A030	3.0m, Male - Male Connectors	
VA15B012	1.2m, Male - Female Connectors	
VA15B015	1.5m, Male - Female Connectors	
VA15B020	2.0m, Male - Female Connectors	
VA15B030	3.0m, Male - Female Connectors	

Gas Sampling Line

- Designed for precision sampling of end-tidal CO₂
- Offered in both male or female Luer for machine fitting
- Different length is available



Closed Wound Drainage System


- 100% medical grade silicone for superior biocompatibility
- Standard system includes evacuator(150ml) and drainage tube
- Available choices: with/without stainless steel trocar, Flat perforated/Flat fluted/Round fluted drainage
- Ensuring efficient drainage and hygienic nursing procedure

Flat Perforated Drainage

Code	Specifications
VD01A0701	7mm x 20cm, flat perforated drainage tube, 150ml evacuator
VD01A1001	10mm x 20cm, flat perforated drainage tube, 150ml evacuator
VD01A0702	7mm x 20cm, flat perforated drainage tube, 150ml evacuator, trocar
VD01A1002	10mm x 20cm, flat perforated drainage tube, 150ml evacuator, trocar



Flat Perforated Type

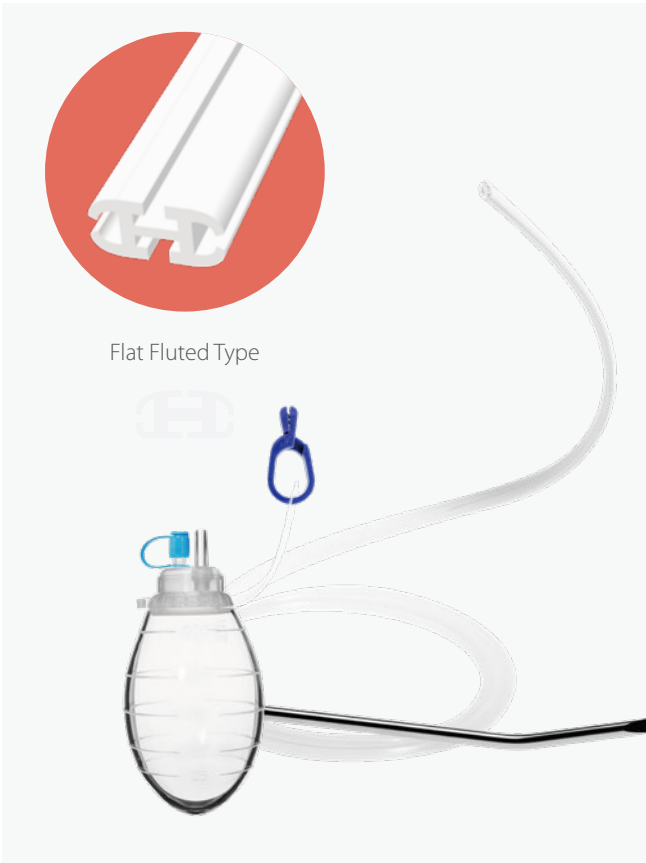

Flat Fluted Type

Flat Fluted Drainage

Code	Specifications
VD02A0701	7mm x 20cm, flat fluted drainage tube, 150ml evacuator
VD02A1001	10mm x 20cm, flat fluted drainage tube, 150ml evacuator
VD02A0702	7mm x 20cm, flat fluted drainage tube, 150ml evacuator, trocar
VD02A1002	10mm x 20cm, flat fluted drainage tube, 150ml evacuator, trocar



Round Fluted Drainage

Code	Specifications
VD03A1001	10Fr x 20cm, round fluted drainage tube, 150ml evacuator
VD03A1501	15Fr x 20cm, round fluted drainage tube, 150ml evacuator
VD03A1002	10Fr x 20cm, round fluted drainage tube, 150ml evacuator, trocar
VD03A1502	15Fr x 20cm, round fluted drainage tube, 150ml evacuator, trocar



Round Fluted Type

Silicone Esophagostomy Feeding Tube

- 100% medical grade silicone for superior biocompatibility
- X-ray opaque line throughout the length
- With 4 side holes
- Markings from the distal end every 5 cm

Code	Specifications	Pack
VG03A0850	8Fr x 500mm	10 pcs/box
VG03A1050	10Fr x 500mm	
VG03A1250	12Fr x 500mm	
VG03A1450	14Fr x 500mm	
VG03A1650	16Fr x 500mm	
VG03A1850	18Fr x 500mm	



TPU Feeding Tube With Guide Wire

- X-ray opaque line throughout the length
- With straight guide wire
- Open distal end with 2 side holes

Code Standard	Specifications	Pack
VG04A0650	6Fr x 500mm	5 pcs/box
VG04A0850	8Fr x 500mm	
VG04A1050	10Fr x 500mm	
VG04A1050	10Fr x 500mm	



Nasogastric Feeding Tube

- X-ray opaque line throughout the length
- With 1cm scaling throughout the tube

Code Standard	Specifications	Pack
VG02A0650	6Fr x 500mm	10 pcs/box
VG02A0850	8Fr x 500mm	
VG02A1050	10Fr x 500mm	



Gastrostomy Feeding Tube

- Gastric balloon for a secure and comfortable fix
- Rounded tip with open distal end
- Rounded suture collar easily adjust to maintain proper tube position

Code Standard	Specifications	Pack
VG05A1201	12Fr, 5ml balloon	1 pc/pk
VG05A1401	14Fr, 5ml balloon	
VG05A1601	16Fr, 5ml balloon	




Extracorporeal Circuit of Hemodialysis

Hemodialysis blood tubing set consists of arterial and venous lines, used during dialysis and are attached with A.V. fistula needle and dialyzer.

- Fluid pathway is smooth inside, no clotting or air bubbles
- Flexible, kink-resistant
- Pump segment with anti-fatigue performance
- Open for tailor-made services

Code	Specifications	Pack
VH05A1201	Blood volume 55ml±10%	10 sets/pk
VH05A1202	Blood volume 26ml±10%	10 sets/pk

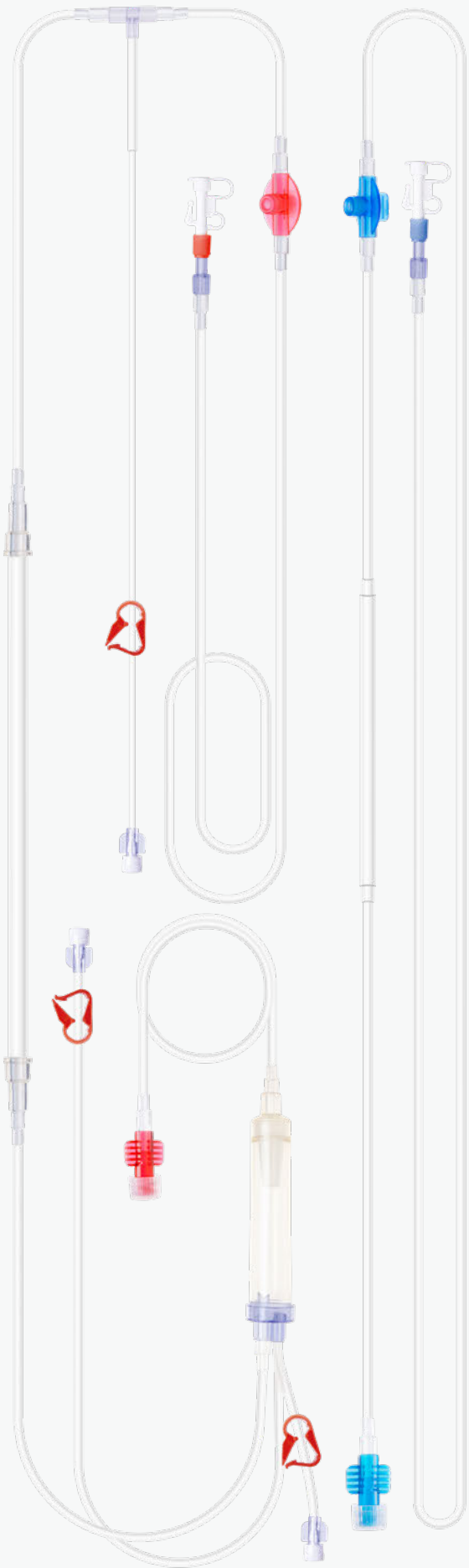

Central Venous Catheter Kit

- Ultra-soft tip prevents risk of vascular trauma
- Made of polyurethane with superior biocompatibility. Long term placement(Up to 28days)
- Radiopacity makes locating more easily after placement
- Scale markings on catheter and guide wire assist in accurate placement
- Standard kit includes: Central Venous Catheter, Steel Guidewire, Dilator, Clamps, Straight-Type introducer Needle with 5ml Syringe, Injection Needle, Injection Caps



Cat. No.	Lumen	Catheter	Guidewire	Dilator	Introducer Needle
A11B011416	Single Lumen 	14Ga x16cm	0.032"x60	8Fr*10cm	18G*6.5cm
A11B011420		14Ga x20cm	0.032"x60	8Fr*10cm	18G*6.5cm
A11B011430		14Ga x30cm	0.032"x70	8Fr*10cm	18G*6.5cm
A11B011616		16Ga x16cm	0.032"x60	6Fr*8cm	18G*6.5cm
A11B011620		16Ga x20cm	0.032"x60	6Fr*8cm	18G*6.5cm
A11B011630		16Ga x30cm	0.032"x70	6Fr*8cm	18G*6.5cm
A11B021813		18Ga x 13cm	0.021"x45	4Fr*5cm	20G*3.8cm
A11B021816		18Ga x 16cm	0.021"x45	4Fr*5cm	20G*3.8cm
A11B021820		18Ga x 20cm	0.021"x60	4Fr*5cm	20G*3.8cm
A11B012013		20Ga x 13cm	0.021"x45	4Fr*5cm	20G*3.8cm
A11B034005	Double Lumen 	4Fr x 5cm	0.018"x45	4Fr*5cm	20G*3.8cm
A11B034008		4Fr x 8cm	0.018"x45	4Fr*5cm	20G*3.8cm
A11B034013		4Fr x 13cm	0.018"x45	4Fr*5cm	20G*3.8cm
A11B034030		4Fr x 30cm	0.018"x65	4Fr*5cm	20G*3.8cm
A11B035008		5Fr x 8cm	0.021"x45	6Fr*8cm	20G*3.8cm
A11B035013		5Fr x 13cm	0.021"x45	6Fr*8cm	20G*3.8cm
A11B035016		5Fr x 16cm	0.021"x45	6Fr*8cm	20G*3.8cm
A11B037016		7Fr x 16cm	0.032"x60	8Fr*10cm	18G*6.5cm
A11B037020		7Fr x 20cm	0.032"x60	8Fr*10cm	18G*6.5cm
A11B037030		7Fr x 30cm	0.032"x70	8Fr*10cm	18G*6.5cm
A11B038011	Triple Lumen 	8Fr x 11cm	0.032"x60	9Fr*10cm	18G*6.5cm
A11B038016		8Fr x 16cm	0.032"x60	9Fr*10cm	18G*6.5cm
A11B038020		8Fr x 20cm	0.032"x60	9Fr*10cm	18G*6.5cm
A11B055508		5.5Fr x 8cm	0.021"x45	6Fr*8cm	20G*3.8cm
A11B055513		5.5Fr x 13cm	0.021"x45	6Fr*8cm	20G*3.8cm
A11B055530		5.5Fr x 30cm	0.021"x60	6Fr*8cm	20G*3.8cm
A11B057016		7Fr x 16cm	0.032"x60	8Fr*10cm	18G*6.5cm
A11B057020		7Fr x 20cm	0.032"x60	8Fr*10cm	18G*6.5cm
A11B057030		7Fr x 30cm	0.032"x70	8Fr*10cm	18G*6.5cm
A11B058516		8.5Fr x 16cm	0.032"x60	9Fr*10cm	18G*6.5cm
A11B058520		8.5Fr x 20cm	0.032"x60	9Fr*10cm	18G*6.5cm
A11B058530		8.5Fr x 30cm	0.032"x70	9Fr*10cm	18G*6.5cm

Mare Uterus Flushing Tube




Intended for use in lavage and cavity washing of the mare's uterus.

Code	Specifications	Pack
VU11A1970	OD 8.7mm, Length 1970mm	5 sets/pk

Obstipation Evacuation Tube

- Inflated cuff at distal end for clearing compacted fecal impaction efficiently
- Clear and transparent tubing with balance of rigidity and flexibility
- Multiple side eyes around the cuff for max irritation
- With atraumatic tip to reduce irritation of inserting
- Needle-free Luer lock hub for better hygienical nursing procedure

Code	Specifications	Pack
VG01A0101	15Fr x 40cm, 3ml Cuff	10 pcs/bag

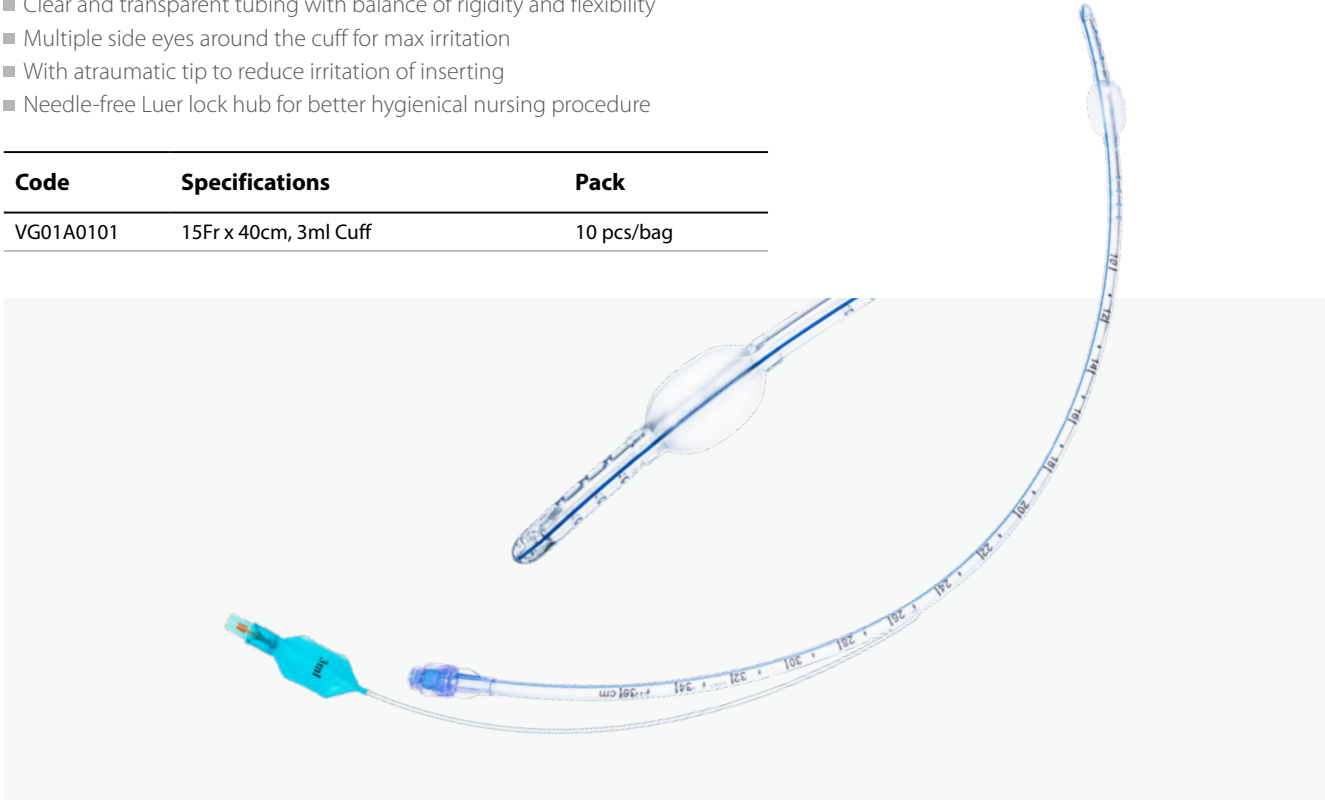

Aerosol Chamber

- Designed for fitting metered dose inhaler
- Efficient medication delivery for pets with asthma, bronchitis, allergies, and breathing difficulties
- With 3 sizes of silicone face masks


Code	Specifications	Pack
VR01A0101	Aerosol chamber set	1 set/box



ETCO₂Sampling Oxygen Cannula

- Designed for oxygen supply and sampling of exhaled CO₂gas simultaneously
- Soft prongs minimize irritation

Code	Specifications	Pack
VR02A0101	Male Luer fitting	10 pcs/bag
VR02B0101	Female Luer fitting	

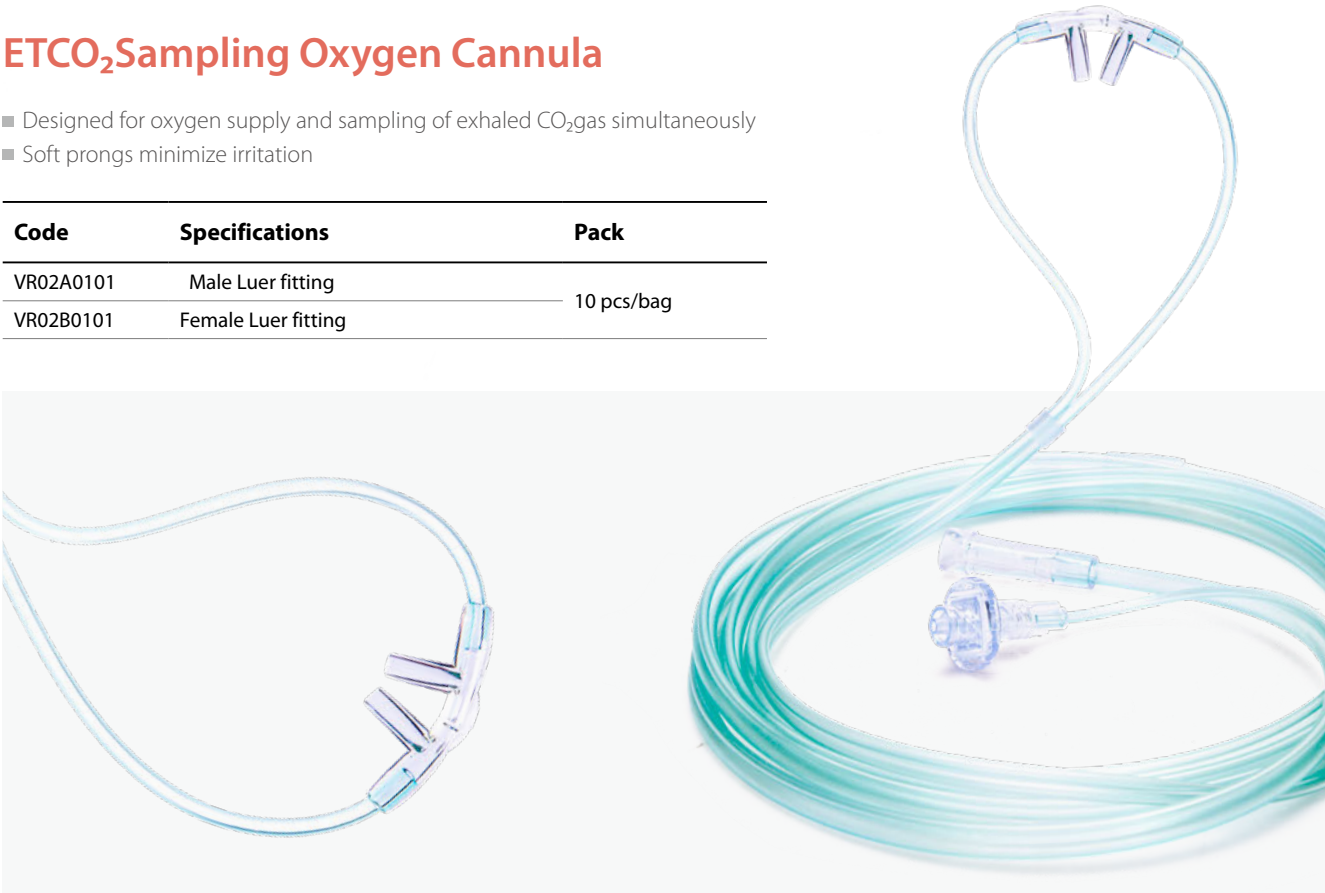

Disposable Yankuer Handle

- The forming angle of handle is designed for more convenient and humanized in clinical use
- Designed for suction needs in various surgical and nursing procedures

Code	Specifications	Pack
VS01A0101	Crown tip	10 pcs/box
VS01B0101	Plain tip	



Suction Connecting Tube with Yankuer Handle

- Transparent material for better visualization
- Striated walls of tube provide superior strength and anti-kinking
- Tubing with great balance of flexibility and rigidity

Code	Specifications	Pack
VS02A0618	Crown tip, 6mm (1/4") x 1.8m (6ft)	10 pcs/box
VS02A0630	Crown tip, 6mm (1/4") x 3.0m (6ft)	
VS02B0618	Plain tip, 6mm (1/4") x 1.8m (6ft)	
VS02B0630	Plain tip, 6mm (1/4") x 3.0m (6ft)	



Suction Tube

- Color-coded connector for easy identification of size
- Soft distal end facilitates comfortable insertion

Code	Specifications	Pack
VS03A0650	6Fr x 500mm	20 pcs/box
VS03A0850	8Fr x 500mm	
VS03A1050	10Fr x 500mm	
VS03A1250	12Fr x 500mm	
VS03A1450	14Fr x 500mm	
VS03A1650	16Fr x 500mm	
VS03A1850	18Fr x 500mm	

