



Pre - Connected Urine Meter

Reduced risk of CAUTI, safeguarding health

Designed to minimize the exposure of the catheterized urinary tract to external contaminants and pathogens, which have been shown to significantly reduce the risk of acquiring Catheter - Associated Urinary Tract Infections (CAUTIs).

Available with latex foley catheter, all silicone foley catheter and hydrophilic foley catheters.

ComboPlus Foley Catheter Tray

Infection control in a completely closed system

The tray offers a streamlined and efficient solution for foley catheterization procedures. It comes equipped with essential supplies to ensure safe and simple insertion while minimizing the risk of cross-infection.

Available for latex foley catheter, all silicone foley catheter and hydrophilic foley catheters.

The standard tray components include:

- A Foley Catheter Preconnected to a Urine Meter
 - Medical Glove
 - Sterile Water Inflation Syringe
 - Lubricating Jelly
 - Wrap
 - Underpad
 - Fenestrated Drape
- Povidon Iodine Swabsticks
 - Plastic Tweezers
 - Alcohol Pad
 - Specimen Container with Label
 - Tray
 - Fixation Device



These components are meticulously selected to meet clinical requirements and can be customized as needed.

Ordering Information




Cat. No.	Description
Urine Meter & Accessories	
F04C022520	1.1m tube, 2.0L bag, hanging straps, needle-free sample port (Luer)
F04C022530	1.5m tube, 2.0L bag, hanging straps, needle-free sample port (Luer)
F04C022540	2.0m tube, 2.0L bag, hanging straps, needle-free sample port (Luer)
F04F012010	2.0L exchange bag with T-tap outlet
F04F013510	3.5L exchange bag with T-tap use-once outlet
F04E010130	Square bed hook
F04E010120	Round bed hook
Pre - Connected Urine Meter	
EB12ADA0001	Silicone foley catheter 12Fr 10mL preconnected to a top-mounted urine meter 2500mL, 1.1m tube, 2.0L bag, hanging straps, needle-free sample port (Luer)
ComboPlus Foley Catheter Tray	
BB12ADA0001	Silicone foley catheter 12Fr 10mL preconnected to a top-mounted urine meter 2500mL, 1.1m tube, 2.0L bag, hanging straps, needle-free sample port (Luer) with standard components

Urine Meter

Diuresis Measuring System

Safe and accurate urine measurement
in critically ill patients

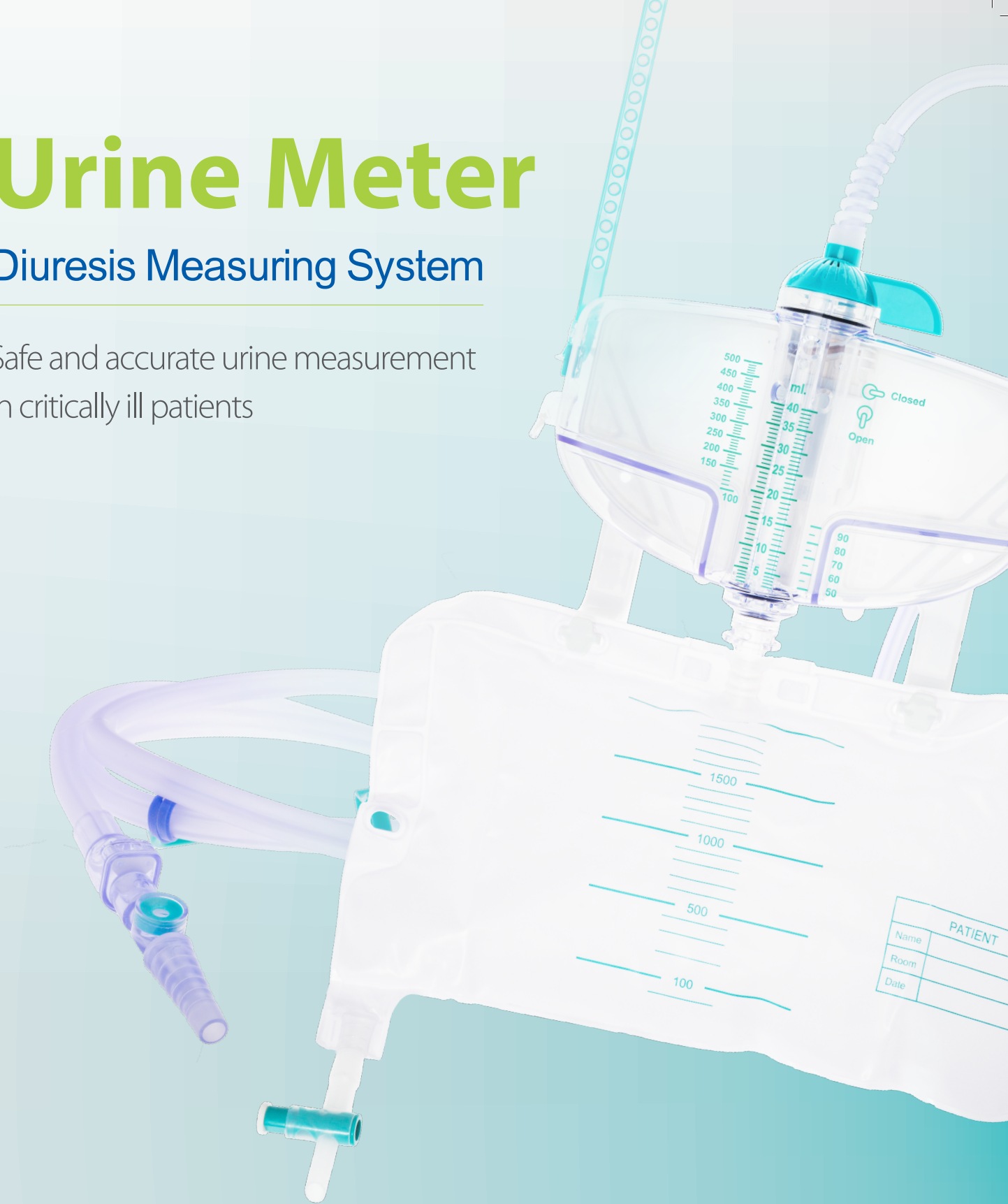

Address: 47 Guomao Avenue South, Panyu, Guangzhou, China
Tel: +86-20-84758878 Fax: +86-20-84758224 Post Code: 511434
Marketing & Trading Office
Add: RM 805, No 55, 21 Long, Pusan Road, Shanghai, China
Tel: +86-21-58397623 Fax: +86-21-58397615 Post code: 200125
E-mail: info@welllead.com.cn www.welllead.com

WELL LEAD MEDICAL CO., LTD.


welllead

Top - Mounted Urine Meter



The urine meter is designed to collect and measure hourly diuresis in critically ill patients by connecting to urinary catheters. The metered chamber precisely measures the urine output, while the drainage bag collects urine, providing indicative volume information.

- 1 Transparent connector with protective cap
- 2 Needle-free sampling port
- 3 1st non-return valve
- 4 Sheep clamp
- 5 Tube sliding clamp
- 6 Double-lumen tubing
- 7 Anti-kink spiral
- 8 Easy-open tap
- 9 Hydrophobic vent filter
- 10 Fine-measuring chamber
- 11 Large chamber with white backing
- 12 Bayonet bag connector
- 13 2nd non-return valve
- 14 Vented collection bag
- 15 Opening
- 16 Large bottom outlet with two stoppers
- 17 2-level hook attachment
- 18 Hanging cord
- 19 Hanging straps

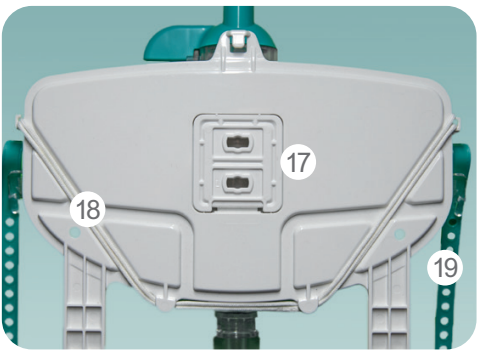

Reliable Measurement

High - precision measurement also with small volume

Fast and Precise Measuring

- The well-defined, dark-colored scaling against the white background allows precise and quick measurement down to 1ml of volume.



Preservation of Fine Measurement

- The measuring chamber design prevents loss of fine measurement when the system is placed on the bed, e.g., during patient transport.



Double - Lumen Tubing with Non - Return Valve

- Splits the urine flow into 2 channels
- Ensures reliable emptying into the chamber without backflow to the patient, also with small volumes.
- Minimizes the risk of urine stasis
- Available Length: 1.1m, 1.5m, 2m



Infection Control

Reliable design against contamination


Top - Mounted Chamber Design

- Minimizes retrograde infection
- Prevents urinary tract infections




2 Non - Return Valves (NRV) Helping to further minimize the risk of retrograde contamination

- The 1st NRV positioned behind the sampling port prevents residual urine from returning in the direction of the bladder.
- The 2nd NRV in the bag prevents fluid reflux back into the chamber.

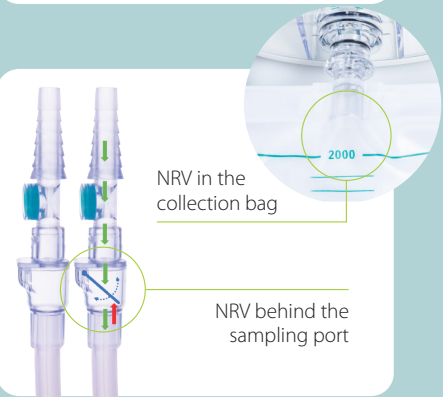

Closed - System with Efficient Venting

- Enhanced with a high-performance bacterial filter for efficient, low-risk venting of the measuring chamber.
- Prevents urine and potential contaminants from exiting the chamber.

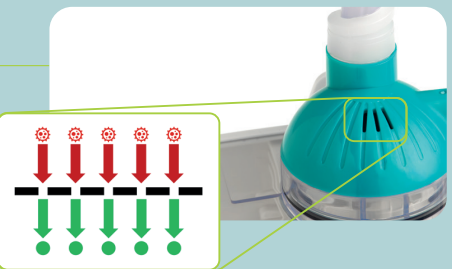

Ease of Use with Safety

One - hand hygienic and practical handling


Easy Empty Tab

- With safe grip allows one - hand operation




Anti - Kink Spiral

- Reduces the risk of blockage and urine stasis

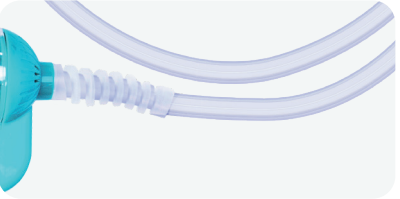

Needle - Free Sampling Port

- Effectively prevents needle - stick injuries

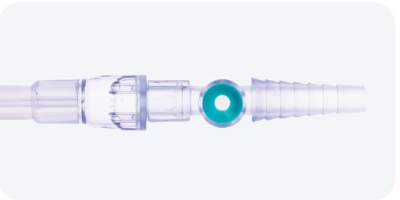

Safety T - Type with Large Outlet

- Ensures one - hand operation with fast empty of the bag



Multiple Hanging Methods For Using on All Types of Beds



Square Hook



Round Hook



Hook



Cord



Straps