

wellead®

WELL LEAD MEDICAL CO., LTD.

Address: 47 Guomao Avenue South, Panyu, Guangzhou, China
Tel: +86-20-84758878 Fax: +86-20-84758224 Post Code: 511434
Marketing & Trading Office
Add: RM 805, No 55, 21 Long, Pusan Road, Shanghai, China
Tel: +86-21-58397623 Fax: +86-21-58397615 Post code: 200125
E-mail: info@wellead.com.cn www.wellead.com



wellead®

ConviCath Intermittent Catheter





ConviCath Intermittent Catheter

**Designed for easy, hygienic and comfortable
intermittent catheterization**

TABLE OF CONTENTS

- 02** What is intermittent catheterization?
- 03** Main advantages at a glance
- 04** Key features of ConviCath intermittent catheters
- 06** Ordering information
- 08** Instruction for use
- 09** Our solution for urinary incontinence

What is intermittent catheterization?


Intermittent catheterization involves inserting a disposable catheter into the bladder to drain the urine. This helps the bladder empty completely. Once the bladder is empty, the catheter is removed and discarded. It may seem intimidating and difficult at first, but it is a fairly simple procedure. Most people are able to catheterize themselves with a little training and practice.

It is recommended intermittent catheterization be performed at regular intervals throughout the day depending on the person's fluid intake and as instructed by the healthcare provider. The ability to perform catheterization and adhere to a schedule is essential to success of the ISC program. You may need to catheterize every 4 - 6 hours to keep the amount of urine in your bladder less than 400-500 ml.

There are many benefits of intermittent catheterization that make it the ideal choice for many individuals, including:

- **Reduce the risk of urinary tract infections (UTIs):** intermittent catheters don't remain in the body like indwelling catheters, so they significantly decrease the risk of infection.
- **More freedom and independence:** it allows you to manage your bladder at times and places that suit you, giving you more control of your life and enhancing your ability to take part in daily activities.
- **Protect bladder function:** emptying your bladder completely prevents residues or urine remaining which will reduce the amount of times you need to visit the toilet. This not only protects the bladder but also reduces the risk of leakage, urinary tract infections and other complications.
- **User friendly:** with practice, many find intermittent catheterization easy. When used properly under the instruction of a healthcare provider, you can safely catheterize on your own.

Main advantages at a glance

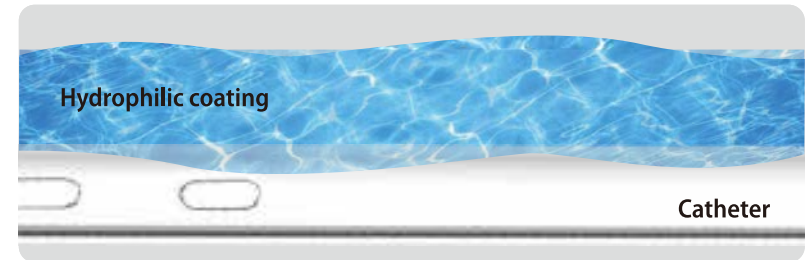

- Hydrophilic coating activated instantly for creating an extremely low friction, slippery surface.
- Ready to use with an integrated water sachet.
- Touch-free hygienic protective catheter sleeve.
- Atraumatic tip and eyes for minimal risk of urethral trauma.
- Different connectors for male and female use. Special ergonomic grip design, tailor made for women's needs.
- Individually packed and sterile, for single use.
- DEHP FREE.



Key features of ConviCath intermittent catheters

Thumb Grip

Allow the package to be easily opened even if manual dexterity is limited. It can also be used to hang the catheter package.

Adhesive Tab

For reliable attachment to any dry and clean surface, allowing both hands to remain free.

Tube

Flexible tube with atraumatic tip and eyes.

Easy-break Water Sachet

With a gentle squeeze to break the sachet. Water comes out in the direction of the arrow to activate the hydrophilic coating

Hydrophilic Coating

Activated easily when the internal water sachet is broken, creating an extremely low friction and slippery surface.

Connector

Different design for male and female use. Special ergonomic grip design is easy to control and is tailor made for women's needs.

Male Connector

Female Connector

Hygienic Protective Catheter Sleeve









Touch free design offers full control with no need to touch the catheter tube during insertion, greatly minimizing the risk of infection.

Ordering information









For Men











Color Code	Size(Fr/Ch)	Pack	Cat.No. w/o water sachet	Cat.No. with water sachet
 Blue	8Fr	box of 30	F01E150810	F01E210810
 Black	10Fr	box of 30	F01E151010	F01E211010
 White	12Fr	box of 30	F01E151210	F01E211210
 Green	14Fr	box of 30	F01E151410	F01E211410
 Orange	16Fr	box of 30	F01E151610	F01E211610
 Red	18Fr	box of 30	F01E151810	F01E211810
 Yellow	20Fr	box of 30	F01E152010	F01E212010
 Violet	22Fr	box of 30	F01E152210	F01E212210




Color Code	Tiemann Tip Size(Fr/Ch)	Pack	Cat.No. w/o water sachet	Cat.No. with water sachet
 White	12Fr	box of 30	F01E171210	F01E241210
 Green	14Fr	box of 30	F01E171410	F01E241410
 Orange	16Fr	box of 30	F01E171610	F01E241610
 Red	18Fr	box of 30	F01E171810	F01E241810
 Yellow	20Fr	box of 30	F01E172010	F01E242010
 Violet	22Fr	box of 30	F01E172210	F01E242210

For Women




Color Code	Size(Fr/Ch)	Pack	Cat.No. w/o water sachet	Cat.No. with water sachet
 Blue	8Fr	box of 30	F01E160810	F01E230810
 Black	10Fr	box of 30	F01E161010	F01E231010
 White	12Fr	box of 30	F01E161210	F01E231210
 Green	14Fr	box of 30	F01E161410	F01E231410
 Orange	16Fr	box of 30	F01E161610	F01E231610
 Red	18Fr	box of 30	F01E161810	F01E231810
 Yellow	20Fr	box of 30	F01E162010	F01E232010
 Violet	22Fr	box of 30	F01E162210	F01E232210

Instruction for use

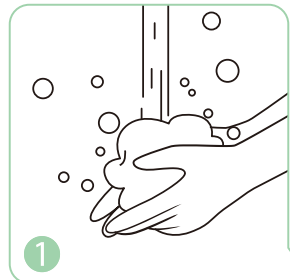


1

Wash hands thoroughly with soap and water.

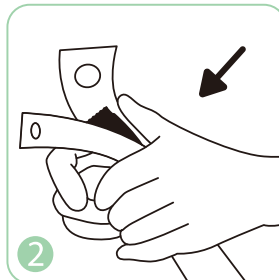


2

Before opening the pouch, squeeze the water sachet and let the water run down into the package.

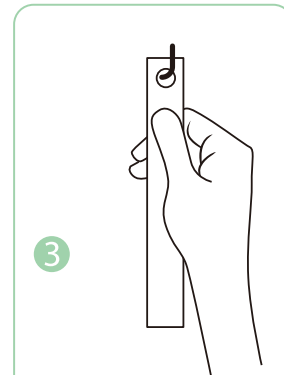


3

Hang the package by its thumb grip or fix it to a dry surface using the adhesive tab.

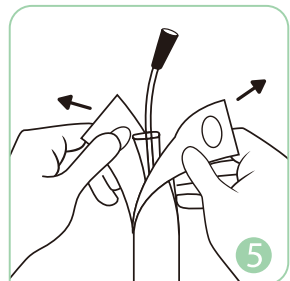


5

Peel the package open by the thumb grips and remove the catheter using the protective sleeve.

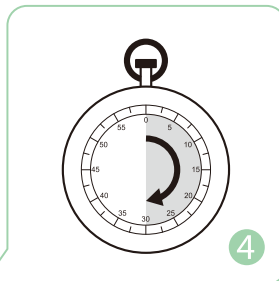


4

Soak the catheter for at least 30 seconds to activate the hydrophilic coating.

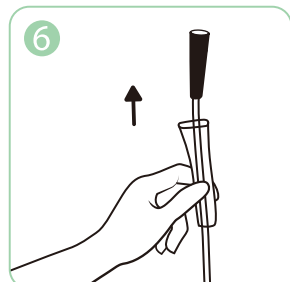

6


Holding the connector, slide the protective sleeve along the catheter towards the tip.

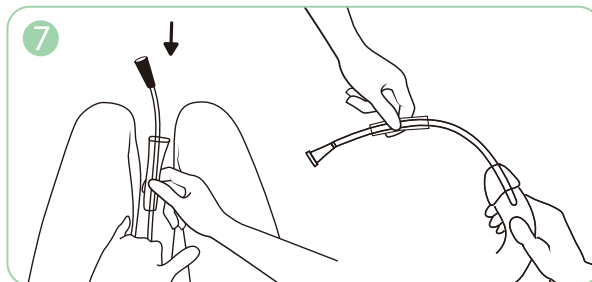

7

Holding the protective catheter sleeve, insert the catheter into the urethra.

Our solution for urinary incontinence

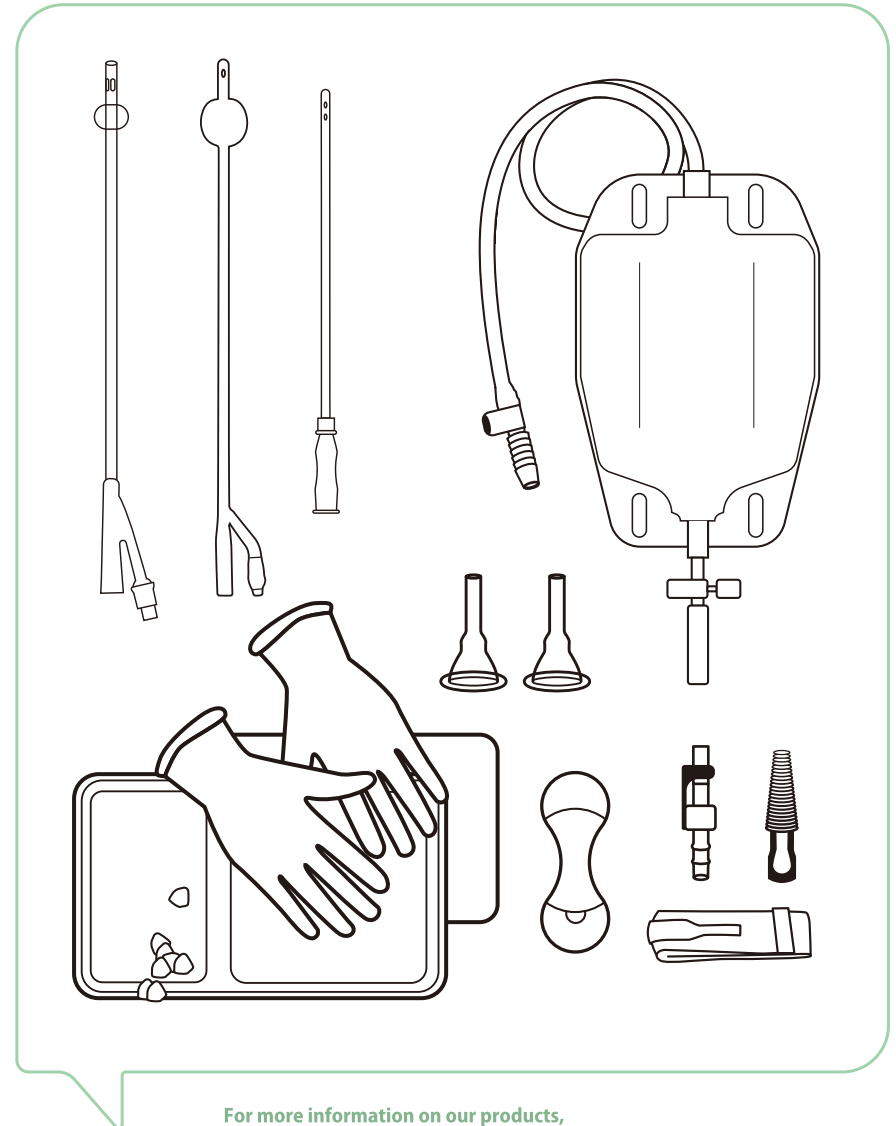

For more information on our products, contact Well Lead Marketing & Trading Office at +86-21-58397623 or visit www.wellead.com