


MADE IN
GERMANY

AUTOCLAVABLE

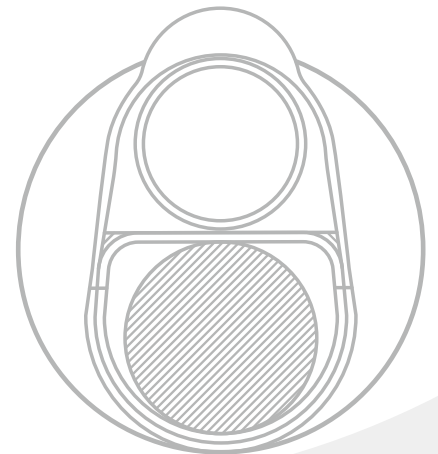

- *Minimum Sheath Size*
- *Maximum Working Channel*

WELL LEAD MEDICAL CO., LTD.

Add: 47 Guomao Avenue South, Panyu Guangzhou, China

Marketing & Trading Office

Add: RM 805, No 55, 21 Long, Pusan Road, Shanghai 200125, China


Tel: +86-021-58397623 Fax: +86-021-58397615

www.welllead.com.cn e-mail: info@welllead.com.cn

wellead®

Features

- Straight, rather than gradual, shaft with minimum diameter for easier, and less invasive advancement;
- Slightly beveled tip for atraumatic insertion;
- Semilunar design providing maximum working channel, and better irrigation flow;
- CPR2/95, CPR2/125, and CPR3/125 allowing the insertion of ureteral Stents;
- Ergonomic hand piece for convenient and comfort control of the device;
- Compatible with all brands of video cameras and light cables;
- Irrigation ports with standard luer lock stopcocks;
- Working channel with standard luer lock ports
- Autoclavable.

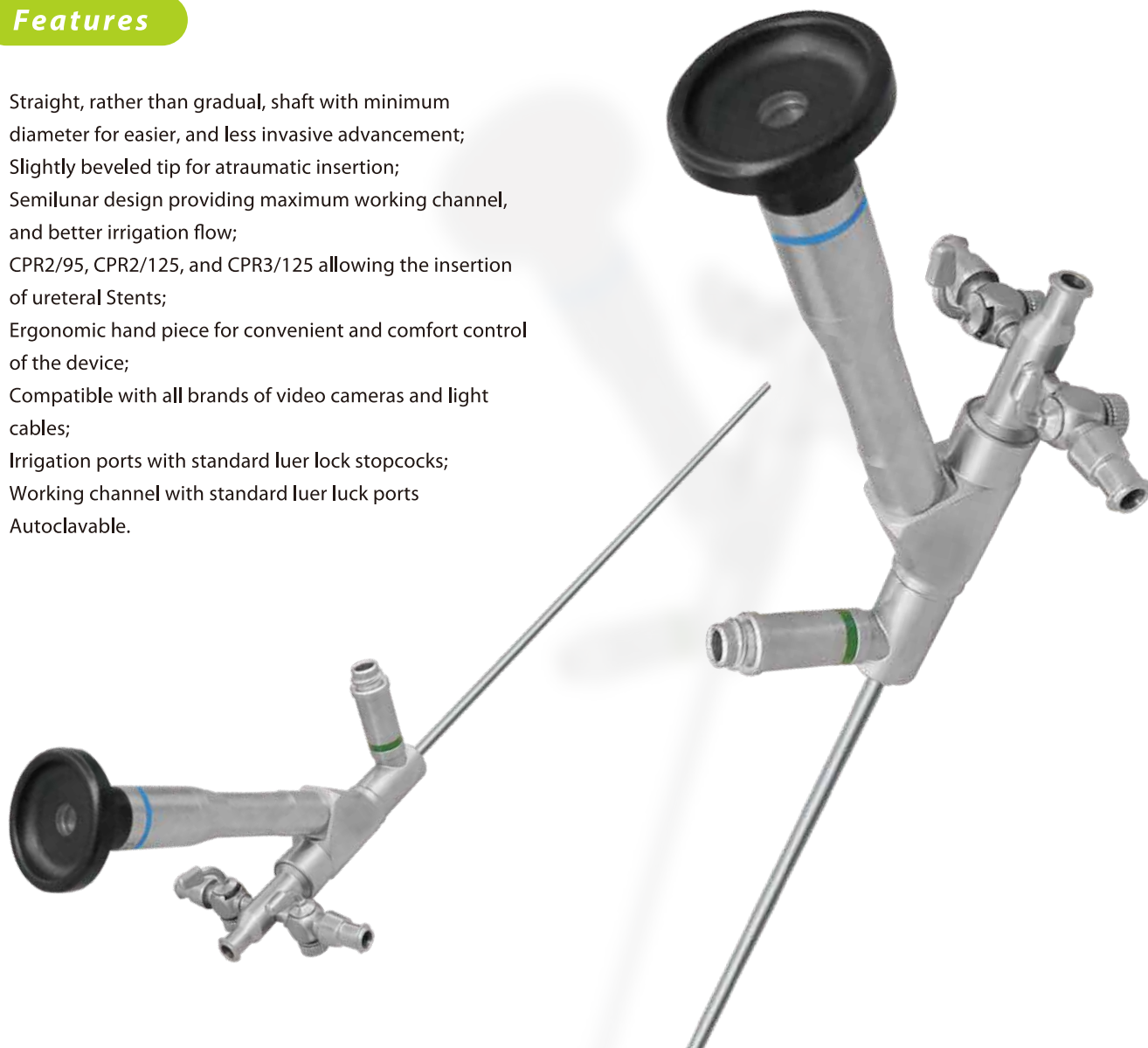

MODEL	URETEROSCOPE	SUPER MINI NEPHROSCOPE	SUPER MINI NEPHROSCOPE	MINI PERC NEPHROSCOPE	CYSTOSCOPE
REF	CPR4/75	CPR2/75	CPR2/95	CPR2/125	CPR3/125
OD	2.5mm (7.5Fr)	2.5mm (7.5Fr)	3.2mm (9.5Fr)	4.2mm (12.5Fr)	4.2mm (12.5Fr)
Working Channel	1.1mm (3.3Fr)	1.1mm (3.3Fr)	1.7mm (5.1Fr)	2.5mm (7.5Fr)	2.5mm (7.5Fr)
Working Length	430mm	250mm	250mm	250mm	340mm
Direction of View	0°	0°	0°	6°	6°
Image	18,000 pixels	18,000 pixels	30,000 pixels	30,000 pixels	30,000 pixels